BEAUTURAL

Percussion Massager

1byone Products Inc.
1230 E Belmont Street, Ontario, CA 91761
Customer Service: +1 909-391-3888
www.1byone.com

Instruction Manual/Bedienungsanleitung/Manuel d'instructions
Manual de instrucciones/Manuale di istruzioni

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before use.

DANGER – To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2) Do not reach for an appliance that has fallen into water. Unplug immediately.
- 3) Do not use while bathing or in a shower.
- 4) Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) This appliance has a polarized plug (one pin is wider than the other). To reduce the risk of electric shock, this plug is intended to fit into a polarized outlet in only one orientation. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

SAVE THESE INSTRUCTIONS.

Thank you for purchasing Beautural Percussion Massager. This instruction manual contains important information about safety, usage, and disposal. Use the product as described and keep this manual for future reference. If you sell this massager or pass it on, also give this manual to the new owner.

Safety Information

This product is only authorized for domestic use and not a substitute for physiotherapeutic or medicinal treatments. It is only intended for applications described in these instructions.

- Under no circumstances may the massager be submerged in fluids, or fluids be permitted to enter the massager's body. Do not expose the massager to high humidity environments. If any liquids enter into the massager's body, disconnect the massager from the power source immediately and allow a qualified technician to inspect it.
- FOR ADULT USE ONLY. This massager is not a toy and should not be used by children. Only allow dependents to use the massager when under proper supervision.
- Do not use this massager before bed. The massage has a stimulating effect and can delay sleep.
- Never use this massager in automobiles.
- Do not operate for more than 15 minutes continuously. It is recommended to turn off the massager and allow a break of 15-20 minutes between each use.
- Remove all scarves, neckties, necklaces and jewelry before using this massager.
- This massager is intended for household use only.

Medical Warnings

The following users should not use this massager or must consult their doctor before use:

- Those who are suffering from osteoporosis or fracture.
- Those suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies.

- Those who are suffering from any of the following: common cold accompanied by fever, varicose veins, thrombosis, phlebitis, jaundice, diabetes, nervous system diseases (e.g. sciatica), tuberculosis, tumors, hemorrhoids, bruises or acute inflammation.
- Those who are suffering from hypertension, heart diseases or brain diseases.
- Those who have a pacemaker or implanted medical device.
- Women who are pregnant or menstruating.
- Those who are under medical treatment or feeling unfit.
- Those who suffering from malignant tumours.
- Those with contagious skin conditions.
- Those with injuries to the body's surface.
- Those who are suffering from acute diseases.
- Those who have recently taken medication or consumed alcohol (resulting in reduced levels of awareness).
- Those who are in bed or sleeping.
- Those under the age of 16.

The massager may not be used:

- on swollen, inflamed or injured body parts.
- for the massaging of bones, especially the spine and the skull.
- close to the eyes or other sensitive body parts.

If you are unsure about your current medical condition, obtain advice from your doctor before using the percussion massager.


Safety Inspection Before Use

Never operate this massager under these circumstances:

- The massager has a damaged cable or plug.
- The massager has been damaged or encountered a heavy fall.
- The massager has been dropped into water.
- The massager generates a burning smell or other odors.
- The power cable or power plug generate abnormal heat.


DO NOT attempt to fix the massager yourself. Contact customer service immediately for further assistance.

Product Description

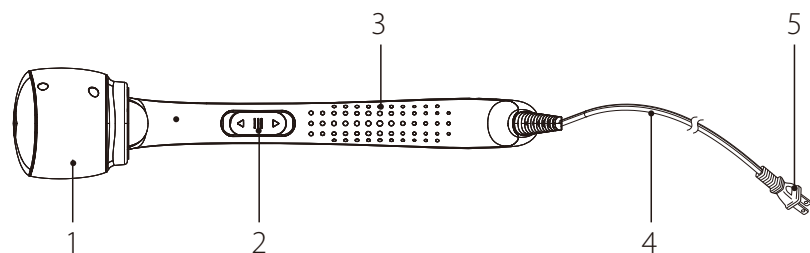- 1. Massaging Node
- 2. OFF/Speed Switch

- 3. Handle
- 4. Power Cable
- 5. Power Plug

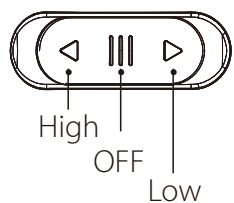Operation

Included Items

- 1 x Percussion Massager
- 1 x Instruction Manual

Before using the massager for the first time, check the package contents to ensure all parts are free of visible damage.

Handling and Operation


1. Always make sure the massager is set to the "OFF" position before inserting the power plug into an electrical outlet.
2. Toggle the OFF/Speed Switch to turn on the massager. Unplug the massager immediately from power after use.
3. Hold the handle and place the massager head lightly on the area you wish to massage.
4. Slowly move the massager head, keeping it over the targeted area. Massage the area for a short period of time before moving to another area. Never massage the same area longer than 3 minutes.

Caution:

- Never use this massager on or near the eyes and other highly sensitive areas.
- Do not apply the massager to the collarbone, pelvis or spine.

Thigh

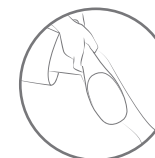Calf

Arm

Waist

Shoulder

Back

- Toggle on OFF/Speed Switch to the OFF position to turn off the massager.
- Toggle on OFF/Speed Switch to the High position for higher intensity.
- Toggle on OFF/Speed Switch to the Low position for lower intensity.

Cleaning & Storage

- Always remove the power plug from the electrical socket before cleaning the massager.
- Unplug the massager and use a dampened soft cloth to clean any dust deposits.
- Never allow water or any other liquid to enter the massager. Use neutral detergent only.
- Protect the massager from light and store it in a cool, dry place.
- Do not hang the massager by the power cable.

Specifications

For US
Rated Voltage: 120V~60Hz
Rated Power: 18W

For EU
Rated Voltage: 230V~50Hz
Rated Power: 18W

Vibration Frequency: 2000-3500RPM
Noise: \leq 55dBA (at 1m/3.3ft)
Size: 430 (Length) x 90 (Width) x 250 (Height) mm
(16.93 x 3.54 x 9.84 inches)
Weight: 1.2Kg (2.6lbs)

Disposal

Disposal of the Appliance

 Under no circumstances should you dispose of the appliance in normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved disposal company or your municipal waste facility. Please observe currently applicable regulations. Please contact your waste disposal center if you need any further information.

 The appliance's packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

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