Rechargeable Neck & Shoulder Shiatsu Massager

1byone Products Inc.
1230 E Belmont Street, Ontario, CA 91761
Customer Service: +1 909-391-3888
www.1byone.com

Instruction Manual/ Bedienungsanleitung
Manuel d'utilisation / Manual de instrucciones / Manuale di istruzioni

IMPORTANT WARNING

- To reduce the risk of electric shock, do not remove the cover. There are no user-serviceable parts inside.
- To reduce the risk of fire or electric shock, do not expose this product to rain or moisture.

The lightning flash with arrow symbol inside an equilateral triangle is intended to alert the user to the presence of uninsulated "hazardous voltage" within the product's enclosure, which may be of sufficient magnitude to induce a risk of electric shock.

The exclamation point inside an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the product.

IMPORTANT: To reduce the risk of electric shock, make sure the plug fits fully into the outlet. If not, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to modify the plug in any way.

IMPORTANT SAFETY INSTRUCTIONS

When using this massager, basic precautions should always be followed.

A. POTENTIAL DANGERS

To reduce the risk of electric shock:


- Always unplug this massager from the electrical outlet immediately after use and before cleaning.
- Do not reach for a massager that has fallen into water. Unplug it immediately.
- Do not use the massager while bathing or showering.
- Do not place or store the massager where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquids.
- Never use pins or other metallic fasteners with this massager.
- Carefully examine the covering before each use. Dispose of the massager if the covering shows any signs of deterioration and causes blistering and racking on the user's skin.
- Keep the massager dry and do not operate it in a wet or moist environment.

B. WARNINGS

To reduce the risk of burns, fire, electric shocks, or injury to persons:

- The massager should never be left unattended when in use.
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shocks, or injury.
- Close supervision is necessary when this massager is used by or near children, dependents, or disabled persons.
- Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this massager under the following circumstances: the massager has a damaged cord or plug; the massager has been damaged or encountered a heavy fall; the massager has been dropped into water.
- Do not carry this massager by its power cord or use the cord as a handle.
- Do not crush or fold this massager when storing away.
- Keep the cord away from heated surfaces.
- Never drop or insert any foreign object into any openings.
- Do not use outdoors. This massager is intended for household and indoor use only.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not operate in the presence of explosive or flammable fumes.
- To disconnect, set all controls to the off position and remove the plug from the outlet.
- Do not overload the electrical outlet. Use only the power outlet type as indicated in the "Specifications" section of this manual.
- To avoid the risk of electric shock, do not disassemble or attempt to repair the massager. Incorrect repair can cause risk of electric shock or injuries when the unit is used.
- Never remove the plug from the outlet by pulling the power cord.

Product Description

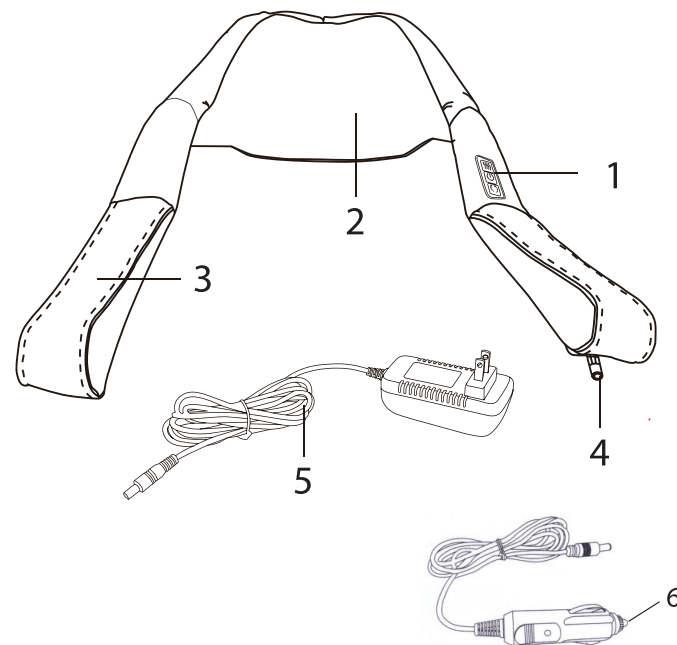- 1. Control Panel
- 2. Shiatsu Nodes
- 3. Positioning Handles

- 4. DC Power Port
- 5. Power Adapter
- 6. Car Adapter

Operating Instructions

1. Check the battery power of the massager.


The massager comes with a built-in rechargeable battery, and can be used immediately upon delivery. A flashing LED on the massager means the unit's battery is running low. In this case, turn off the massager and charge using the enclosed power adapter. The flashing LED will become solid again when the massager is fully charged.

Note: The massager can be used during charging.

2. Grip the positioning handles and drape the main unit around your neck with the Shiatsu nodes at the center of your back.


3. Press and hold the  ON/OFF Button for 3 seconds and the massager immediately starts in Auto Reversing Mode. Default Shiatsu intensity will be high.

4. Press the  Intensity Setting Button to adjust the intensity. The default intensity will be high. Press the Intensity Setting Button to advance the massager to Medium intensity, and press it again for Low intensity.
5. Press the  Auto Reversing Mode Button to switch the Shiatsu Nodes between Auto Reversing Mode and Single Direction Mode. The default mode is Auto Reversing Mode, which will change mode direction every minute. every one minute.
6. Press the  Heat Button to activate heat. Press it again to deactivate heat.
7. Use the positioning handles to reposition the Shiatsu nodes on top of the needed treatment area.
8. Press the  ON/OFF button at any time to turn off the massager.
9. The massager will shut off automatically after 15 minutes of continuous use.

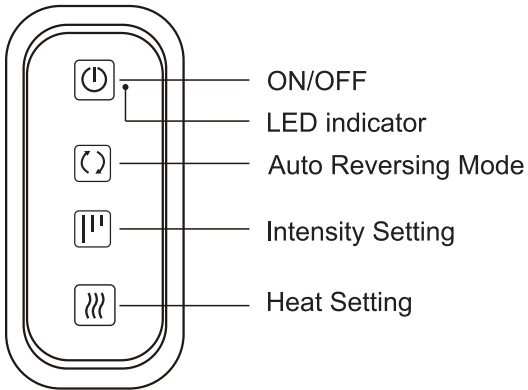Note: It is NOT recommended to use the heat function for more than 15 minutes at the time. Allow the massager to cool down between heat cycles.

Treatment Areas

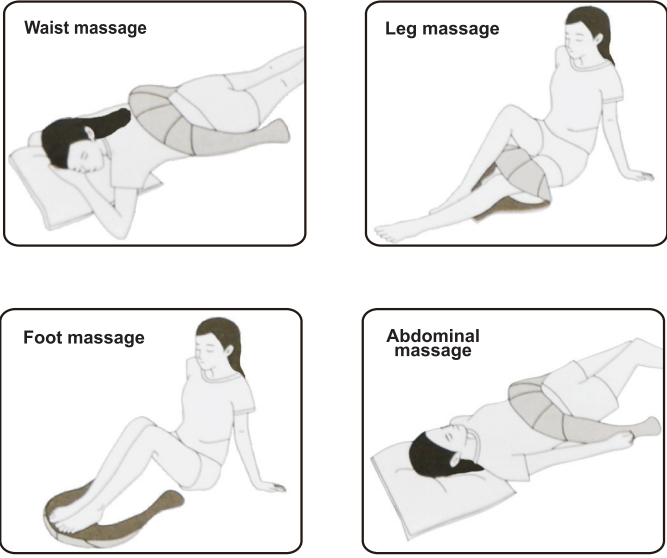Maintenance and Care

1. Place the massager in a safe, cool, and dry place when not in use.
2. After unplugging, wipe with a soft, damp cloth using non-abrasive cleaners.
3. Never submerge the massager in any liquids.
4. Keep away from all solvents and harsh detergents.
5. DO NOT attempt to repair this massager. There are no user-serviceable parts.

Specifications

Input: 100V-240V~ 50/60Hz
Output: 14V 1.7A
Battery Capacity: 11.1V, 2200mAh
Charging Time: Approximately 3 hours
Continuous Use Time: 1~1.5 hours

Disposal

Disposal of the Appliance

Under no circumstances should you dispose of the appliance in normal domestic waste.

This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved disposal company or your municipal waste facility. Please observe currently applicable regulations. Please contact your waste disposal center if you need any further information.

Disposal of Batteries

Used batteries may not be disposed of in household waste.

All consumers are statutorily required to dispose of batteries at a collection point provided by their local municipality or retail store. The purpose of this obligation is to ensure batteries are disposed of in a non-polluting manner. Only dispose of batteries when they are fully discharged.

The appliance's packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

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